\$FLIR



mobile GC-MS with real-time air monitoring FLIR GRIFFIN $^{\!\scriptscriptstyle \mathrm{TM}}$ G465

The FLIR Griffin G465 is the most flexible mobile GC-MS (Gas Chromatograph/ Mass Spectrometer) system available for mobile and forward deployed labs, reconnaissance vehicles, and other portable platforms. It is equipped with a rugged, internal shock isolation system that is tested to rigorous MIL-STD-810G standards. The Griffin G465 contains the same standard injection port commonly found on lab-based GC-MS systems and accepts sample injections via PSI-Probe™, syringe, SPME fiber, headspace, and autosampler. It also includes an additional integrated universal sampling port, providing direct air sampling capability, while also accepting plug-and-play samplers like the Griffin X-Sorber and Purge and Trap. The addition of a heated sampling line enables survey mode for near real-time vapor threat identification. The Griffin G465 accurately detects and identifies explosives, drugs, CWAs, TICs, environmental pollutants, and other chemicals. The simplified user interface gives field-operators and scientists quick and accurate answers. The Griffin G465 helps preserve sample integrity, eliminate the expense of shipping unwanted samples, and lead to real-time, actionable countermeasures that protect public safety.

www.flir.com/G465



SURVEY MODE FOR RAPID VAPOR THREAT ID

Vehicle-ready GC/MS platform for mission-critical rapid chemical identification

- Membrane-based, near real-time vapor threat detection with integrated heated sampling line
- Internal active vacuum pumping maximizes system uptime and expands range of sampling options
- Integrated shock isolation system minimizes footprint on a vehicle
- MIL-STD-810G compliant for shock, vibration, temperature, and humidity
- Optional vehicle mounting kit simplifies installation for operation in a moving vehicle



MULTI-MISSION SAMPLING

Integrated inlets and interchangeable sampling tools provide in-field flexibility and long-term value

- Flexible sampling options accommodate solid, liquid, and vapor samples
- Same split/splitless injector found on typical lab-based GC-MS systems and accepts PSI-Probe™, direct syringe, autosampler, SPME, and headspace injections
- Built-in air sampling module provides 24/7 continuous air monitoring capability while providing a docking station for portable plugand-play water and vapor samplers
- Versatile, quick-connect sampling tools extend the mission to remote locations and eliminate sample prep at the site of action



CONFIDENT, ACTIONABLE CHEMICAL IDENTIFICATION

Industry-leading chemical searching techniques present a single identification result

- Easy-to-use wizard tool guides user through method selection
- Go/no-go alarms eliminate data interpretation
- Chemical ID in five to fifteen minutes
- Access to detailed chromatograph and mass spectral data
- Ability to create mission-specific methods and chemical libraries
- · Operator training completed in one day

SPECIFICATIONS

Griffin G465

Technology GC/MS; fully integrated low thermal mass gas chromatograph (LTM-GC) and MS/MS-capable ion trap mass analyzer

Mass Range / Scan Rate 35-425 m/z; up to 10,000 m/z per second @ 20 points per m/z

Ionization Type

Internal electron ionization (EI) Conversion dynode electron multiplier

ITM-GC Column

Detector

Standard VB-5MS (15M x .18mm x .18µm); others available;

programmable 40 to 300 °C

Calibrant

On-board FC-43 (Perfluorotributylamine)

Carrier Gas

Connection for external gas source (choice of He or H2); gas available from many vendors; H2 generator available

Sampling & Analysis

Sample Introduction

Split/splitless injector accepts:

- Direct syringe injection (1 syringe included)
- SPME fiber (optional)
- Manual headspace sampler (optional)
- Autosampler (optional)
- PSI-Probe[™] thermal separation via TAG[™] (optional)
- PSI-Probe thermal separation via GERSTEL-Twister® (optional)

Universal Sampling Port with standard Tenax TA and Carboxen 1017 dual bed preconcentration tubes accepts:

- Direct air intake onto precon tubes (adaptor included)
- Griffin Purge and Trap for water analysis (optional)
- Griffin X-Sorber handheld vapor sampler (optional)

Membrane Introduction Mass Spectrometry (MIMS) Inlet:

- Direct air intake for near real-time air monitoring capability

Sample Phase Solid, liquid, and vapor

Threats

Detects and identifies explosives, narcotics, CWAs, TICs,

environmental pollutants, and other chemicals

Sampling & Analysis Full identification in 5-15 minutes for most chemicals; near

real-time in survey mode

System Interface

Display & Alerts Full automation by connection with computer

Communication

Ethernet connection TCP/IP; remote operation and diagnostics

Data Storage

Data automatically stored on supplied laptop (500 gb)

Simplified User Interface

Griffin System Software (GSS); GriffinLib, NIST, and AMDIS mass spectral libraries included; capable of user-defined library

Training Requirements

1-2 days depending on level of training desired; Operator, Developer, and Full System certifications available

Specifications are subject to change without notice. For the most up-to-date specs, go to www.flir.com

Power

100-240 VAC; 24 VDC (+/- 5%, 25 A, 1000 W) Input Voltage Cold Start Time <30 mins (includes automatic tuning/calibration)

-13 to 131 °F (-25 to 55 °C)

Environmental

Operating Temp / Humidity

41 to 104 °F (5 to 40 °C); <85% relative humidity

Storage Temp

Physical Features Dimensions (L x W x H)

19.2 x 19.2 x 21.1 in (48.8 x 48.8 x 53.6 cm)


Weight

Enclosure & Protection

99.5 lbs (45.0 kg)

Rugged, internal shock mounting system; integrated vacuum system contains mini turbomolecular pump and quad diaphragm;

no external shock table or vacuum system required

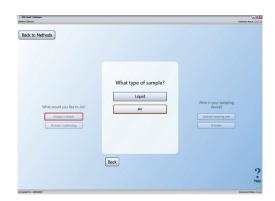



HEADQUARTERS

FLIR Systems, Inc. 27700 SW Parkway Ave. Wilsonville, OR 97070 USA

DETECTION SALES, **AMERICAS**

FLIR Systems, Inc. 1201 S. Joyce Street Suite C006 Arlington, VA 22202 PH: +1-877-692-2120

DETECTION SALES, APAC

FLIR Detection, Inc. 10 Kallang Avenue #09-10 Aperia Tower 2 Singapore 335910 PH: +65-6822-1596

DETECTION SALES, EMEA

FLIR Detection, Inc. Luxemburgstraat 2 2321 Meer Belgium PH: +32 (0) 3665 5106 detection@flir.com

www.flir.com NASDAQ: FLIR

Equipment described herein is subject to US export regulations and may require a license prior to export. Diversion contrary to US law is prohibited. Imagery for illustration purposes only. Specifications are subject to change without notice. ©2020 FLIR Systems, Inc. All rights reserved. 05/20/20

20-0749-DET-DET-DATASHEET-REV-G465 LTR



The World's Sixth Sense®